

Catalogue No.

 AB0037-200
 AB0037-500

Qty:

 400 µg
 1 mg

Anti-CANX

Source: Goat

General description: Goat polyclonal to CANX (Calnexin) - endoplasmic reticulum (ER) membrane marker. CANX is a member of the Calnexin family of molecular chaperones. This protein is a calcium-binding, ER-associated protein that interacts transiently with newly synthesized N-linked glycoproteins, facilitating protein folding and assembly. It may also play a central role in the quality control of protein folding by retaining incorrectly folded protein subunits within the ER for degradation.

Alternative names: Calnexin, CALX, CNX, FLJ26570, histocompatibility complex class I antigen binding protein p88, IP90, major histocompatibility complex class I antigen-binding protein p88, MS952, P90 antibody.

Form: Polyclonal antibody supplied as a 200 or 500 µl (2 mg/ml) aliquot in PBS, 20% glycerol and 0.05% sodium azide. This antibody is epitope-affinity purified from goat antiserum.

Immunogen: Purified recombinant peptide within residues 550 aa to the C-terminus of human CANX produced in E. coli.

Specificity: Detects a band of 90 kDa by Western blot in the following human (293A, primary fibroblasts, HaCat, HeLa, HMEC-1, Jurkat, MNT1, U-118, rat (TR-iBRB), mouse (3T3, AtT-20, Hepa, Raw264.7), monkey (COS-7) and canine (D17) whole cell lysates.

Reactivity: Reacts with Human, Rat, Mouse, Monkey and Canine proteins

Sample	WB	IHC (F)	IHC (P)	IF	ELISA
Human	+++	+++	+++	+++	ND
Rat	+++	+++	+++	+++	ND
Mouse	+++	+++	+++	+++	ND
Canine	+++	+++	+++	+++	ND
Monkey	+++	+++	+++	+++	ND

Usage:


WB: 1:500-1:5,000
 IHC (F): 1:200-1:1,000
 IHC (P): 1:200-1:1,000
 IF: 1:50-1:500

Storage: For continuous use, store at 2-8 C for one-two days. For extended storage, store in -20 C freezer. Working dilution samples should be discarded if not used within 12 hours.


Special instructions: The antibody solution should be gently mixed before use. Avoid freeze/thaw cycles..

References:


1. Saludas L, Garbayo E, Ruiz-Villalba A, et al. *Eur J Pharm Biopharm* 2022 Jan PMID: 34968647
2. Schulte A, Bieniussa L, Gupta R, et al. *Cell Calcium* 2022 Jan PMID: 34896701
3. Saviano A, De Vita S, Chini MG, et al. *Biomolecules* 2021 Oct PMID: 34680124
4. Saludas L, Garbayo E, Ruiz-Villalba A, et al. *Eur J Pharm Biopharm* 2021 Dec PMID: 34968647
5. Ghanawi H, Hennlein L, Zare A, et al. *Nucleic Acids Res* 2021 Dec PMID: 34850154
6. Monteiro-Alfredo T, Caramelo B, Arbeláez D, et al. *Nutrients* 2021 PMID: 34578832
7. Monteiro-Alfredo T, Oliveira S, Amaro A, et al. *Nutrients* 2021 Aug PMID: 34445015
8. Selezneva A, Yoshida M, Gibb A, et al. *Pharmacological Reports* 2021 PMID: 33713314
9. Martinho JM, Coelho A, Oliveiros B, et al. *Int Endod J* 2021 May PMID: 33999433
10. Tavares G, Marques D, Barra C, et al. *Mol Metab* 2021 Feb PMID: 33933677
11. Garbayo E, Ruiz-Villalba A, Hernandez SC, et al. *Acta Biomaterialia* 2021
12. Rodrigues T, Borges P, Mar L, et al. *Pharmacol Res* 2020 Sep PMID: 32942016
13. Sousa FJ, Correia RG, Cruz AF, et al. *Brain, Behavior, & Immunity - Health* 2020 May
14. Monteiro-Alfredo T, Matafome P, Iacia BP, et al. *Oxid Med Cell Longev* 2020 Mar PMID: 32256951
15. Tejero R, Balk S, Franco-Espin J, et al. *iScience*. 2020 Feb PMID: 31981925
16. Graça SSF. MSc Thesis, Universidade do Algarve, Portugal 2019
17. Balk SM, PhD Thesis University of Wurzburg, Germany 2019
18. Fricke S, Metzendorf K, Ohm M, et al. *Cell Rep* 2019 Oct PMID: 31618635
19. Neves C, Rodrigues T, Sereno J, et al. *Oxid Med Cell Longev* 2019 Jun PMID: 31341532
20. Fonseca LMO, MSc Thesis, University of Coimbra, Portugal 2018
21. Silva MM, Gomes-Alves P, Rosa S, et al. *J Biotechnol*. 2018 Aug PMID: 30165116
22. Rodrigues TDA, PhD Thesis, University of Coimbra, Portugal 2018
23. Ribeiro M, Castelhana J, Petrella LI, et al. *J Magn Reson Imaging* 2018 Jan PMID: 29377412
24. Awadh A, PhD Thesis, University of Alberta, Canada 2018
25. Ribeiro STF, PhD Thesis University of Lisbon, Portugal 2017
26. Cabral AMD, MSc Thesis, University of Lisbon 2017
27. Rodrigues T, Matafome P, Sereno J, et al. *Sci Rep* 2017 May PMID: 28490763
28. Thieleke-Matos C, Lopes da Silva M, Cabrita-Santos L et al. *Cellular Microbiology* 2016 Mar PMID: 26399761
29. Sivadasan R, PhD Thesis, Wurzburg University 2016
30. Neves CAF, MSc Thesis University of Aveiro, Portugal 2015
31. Alves LS, Marques ARA, Padrao N, et al. *J Cell Sci* 2021 Sep. PMID: 34528688
32. Gil-Cabrerizo P, Saludas L, Felipe Prósper F et al. *Int J Pharm* 2022 Dec PMID: 36332831
33. Tavares G, Rosendo-Silva D, Simões F, et al. *Int J Mol Sci* 2023 Jan PMID: 36768789
34. Amaro A, Sousa D, Sá-Rocha M, et al. *Life Sci* 2023 Mar PMID: 36948389

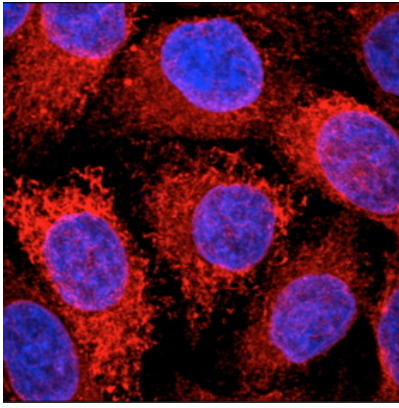


Immunofluorescence – anti-CANX Ab in HeLa cells at 1/250 dilution; cells were fixed with 4% of PFA;

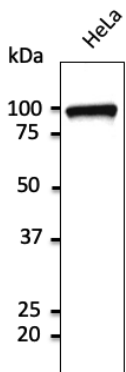



Anti-Calnexin - ER membrane marker Ab at 1/2,500 dilution; lysates at 50 μ g per lane; Rabbit polyclonal to goat IgG (HRP) at 1/10,000 dilution;

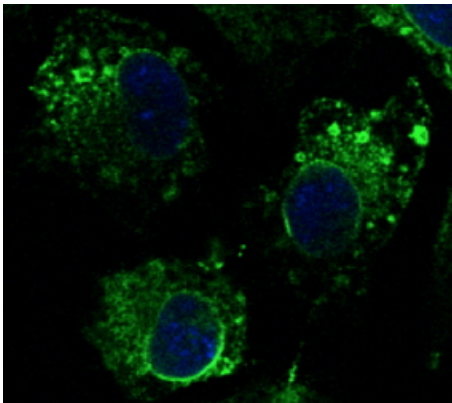



Immunofluorescence – anti-CANX Ab in Hepa1-6 cells at 1/100 dilution; cells were fixed with 4% of PFA;

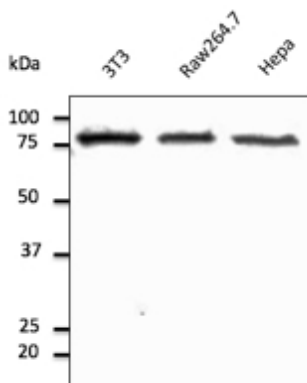



Anti-Calnexin - ER membrane marker Ab at 1/2,500 dilution; lysates at 50 μ g per lane; Rabbit polyclonal to goat IgG (HRP) at 1/10,000 dilution;

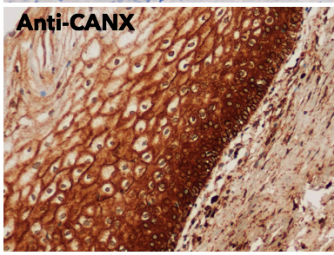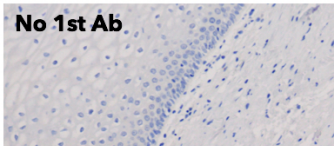



IHC of human cervix using CANX antibody and FFPE tissue after heat-induced antigen retrieval. Anti-CANX Ab at 1:500/DAB detection.

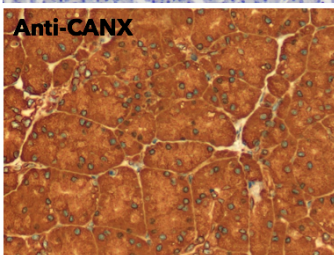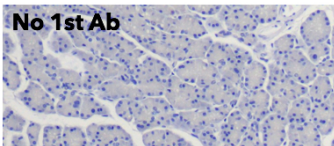



IHC of human pancreas using anti-CANX antibody and FFPE tissue after heat-induced antigen retrieval. Anti-CANX Ab at 1:500/DAB detection.

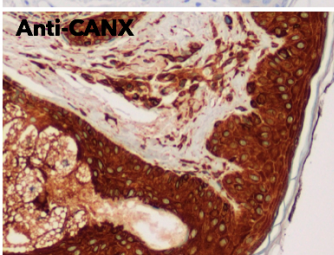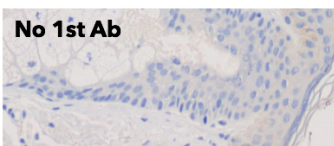



IHC of human skin using anti-CANX antibody and FFPE tissue after heat-induced antigen retrieval. Anti-CANX Ab at 1:500/DAB detection.

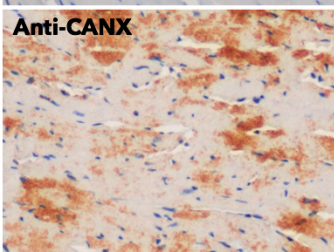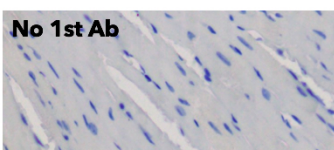

IHC of mouse myocardium using CANX antibody and FFPE tissue after heat-induced antigen retrieval. Anti-CANX Ab at 1:500/DAB detection.

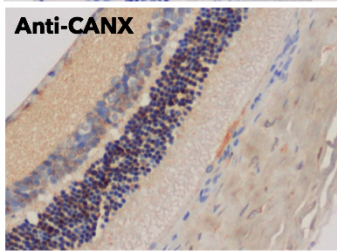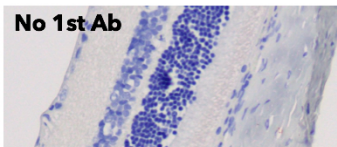

IHC of mouse eye using anti-CANX antibody and FFPE tissue after heat-induced antigen retrieval. Anti-CANX Ab at 1:500/DAB detection.

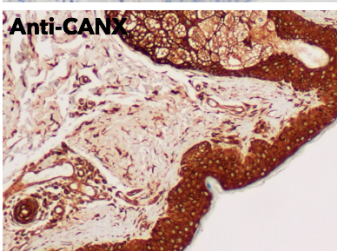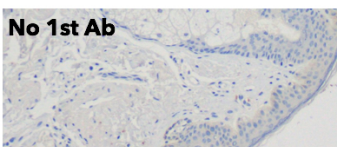

IHC of human skin using anti-CANX antibody and FFPE tissue after heat-induced antigen retrieval. Anti-CANX Ab at 1:500/DAB detection.

For research use only, not for diagnostic use

SICGEN's Proprietary Immunogen Policy

In order to produce high specific antibodies SICGEN has invested a lot of time and effort into selecting immunogen sequences. SICGEN has decided to protect this information by not publishing it on the website. However, these sequences are available on request.