

**Catalogue No.****Qty:**

AB0020-200

400 µg

AB0020-500

1 mg

**Anti-GFP****Source:** Goat

**General description:** Goat polyclonal antibody to GFP (green fluorescent protein). GFP is a protein composed of 238 amino acid residues (26.9 kDa) that exhibits bright green fluorescence when exposed to blue light. In cell and molecular biology, the GFP protein is frequently used as a reporter of expression.

**Alternative names:** Green fluorescent protein antibody.

**Form:** Polyclonal antibody supplied as a 200 or 500 µl (2 mg/ml) aliquot in PBS, 20% glycerol and 0.05% sodium azide. This antibody is epitope-affinity purified from goat antiserum.

**Immunogen:** Purified recombinant peptide produced in E. coli.

**Specificity:** In 293HEK cells transfected with cds plasmid detects a band of 27 kDa by Western blot. This antibody does recognize YFP and CFP and does not recognize RFP and mCherry fluorescent proteins.

**Reactivity:** Reacts with Transfected cells proteins

Sample	WB	IHC (F)	IHC (P)	IF	ELISA
Transfected cells	+++	+++	+++	+++	ND

+++ excellent, ++ good, + poor, ND not determined

**Usage:**

WB: 1:500-1:5,000

IF: 1:50-1:1,000

IHC (P): 1:50-1:1,000


IHC (F): 1:50-1:1,000

**Storage:** Store at -20 C for long-term storage. Store at 2-8 C for up to one month.



**Special instructions:** The antibody solution should be gently mixed before use. Avoid freeze/thaw cycles..

**References:**


1. Cabaço LC, Bento-Lopes L, Neto MV, et al. *JID Innov* 2022 Jun PMID: 36090299
2. Yilmaz V, Louca P, Potamiti L, et al. *Elife* 2022 May PMID: 35533001
3. Martins-Marques T, Costa MC, Catarino S, et al. *EMBO Rep* 2022 May PMID: 35593040
4. Ferreira JV, Soares AR, Ramalho J, et al. *Sci Adv* 2022 Mar PMID: 35333565
5. Wang WA , Demaurex N, et al. *J Biol Chem* 2022 Mar PMID: 35065962
6. Escrevente C, Bento-Lopes L, Ramalho JS, et al. *J Cell Sci* 2021 May PMID: 34002205
7. Annibaldis G, Domanski M, Dreos R, et al. *Nucleic Acids Res* 2020 Sep PMID: 32941650
8. Shoaito H, Chauveau S, Gosseaume C, et al. *J Cell Mol Med* 2020 Jun PMID: 32519441
9. Panayiotou T, Michael S, Zaravinos A, et al. *PLOS Pathog* 2020 Apr PMID: 32298395
10. Wu JS, Yi E, Manca M, et al. *Elife* 2020 Jan PMID: 31975688
11. Alvarez C, Campos AC, Howard JC, et al. *Methods Mol Biol.* 2020 PMID: 31758463
12. Carey LM, PhD Thesis, Indiana University, US 2019
13. Zhang L, Noguchi YT, Nakayama H, et al. *Cell Rep* 2019 Nov PMID: 31747590
14. Fukuda S, Kaneshige A, Kaji T, et al. *Elife* 2019 Sep PMID: 31545169
15. Matsuda T, Oinuma I. *Sci Rep* 2019 Aug PMID: 31383899
16. Geller S, Arribat Y, Netzahualcoyotzi C, et al. *Cell Metab* 2019 Aug PMID: 31474567
17. Zhang L, Uezumi A, Kaji T, et al. *Int J Mol Sci* 2019 Jul. PMID: 31277245
18. Alenquer M, Vale-Costa S, Etibor TA, et al. *Nat Commun* 2019 Apr PMID: 30967547
19. Koike M, Yutoku Y, Koike A. *FEBS Open Bio* 2019 Jan PMID: 31115163
20. Vyas P, Wu JS, Jimenez A, et al. *Sci Rep* 2019 Apr PMID: 30944354
21. Alenquer M, Vale-Costa S, Sousa AL, et al. *bioRxiv* 2018
22. Larson VA, Mironova Y, Vanderpool KG, et al. *Elife* 2018 Mar PMID: 29596047
23. Festa BP, Chen Z, Berquez M, et al. *Nat Commun* 2018 Jan PMID: 29323117
24. Ramos-Nascimento A, Kellen B, Ferreira F, et al. *J Cell Sci* 2017 Dec PMID: 29061883
25. Wu JS, Vyas P, Glowatzki E, et al. *J Comp Neurol* 2017 Oct PMID: 29055051
26. Agarwal A, Wu PH, Hughes EG, et al. *Neuron* 2017 Feb PMID: 28132831
27. Mimura Y, Takagi M, Clever M, et al. *J Cell Sci* 2016 Nov PMID: 27802161
28. Vyas P, Wu JS, Zimmerman A, et al. *J Assoc Res Otolaryngol* 2016 Sep PMID: 27696081
29. Vale-Costa S, Alenquer M, Sousa AL, et al. *J Cell Sci* 2016 Mar PMID: 26940915
30. Martins-Marques T, Anjo SI, Pereira P, et al. *Mol Cell Proteomic* 2015 Aug PMID: 26316108
31. Sargiannidou I, Kim GH, Kyriakoudi S, et al. *Neurogenetics* 2015 Jul PMID: 25771809
32. Kourti M, Ikonomou G, Giakoumakis NN, et al. *Cell Signal* 2015 Mar PMID: 25744540
33. Monteiro NCO, MSc Thesis, University of Coimbra, Portugal 2015
34. Ferreira RRS, MSc Thesis, University of Coimbra, Portugal 2015
35. Paukert M, Agarwal A, Cha J, et al. *Neuron* 2014 Jun PMID: 24945771
36. Moreiras HAF, MSc Thesis, University of Lisbon, Portugal 2014
37. Pais SV, Milho C, Almeida F, et al. *PLoS One* 2013 Feb PMID: 23431368
38. Matlac DM, PhD Thesis, University of Lübeck, Germany 2021
39. Watson EO, MRes Thesis, Bristol University, UK 2022
40. Mirfakhar F, Castanheira J, Domingues R, et al. *bioRxiv* 2023
41. Veltheim M, Master Thesis, Tampere University, Finland 2023
42. Yang Y, PhD Thesis, Tubingen , Germany2023
43. Duarte F, Ramosaj M, Hasanovic Ed, et al. *Mol Ther Methods Clin Dev* 2023 Mar PMID: 37090478
44. van Bruggen R, Patel ZH, Wang M, et al. *eNeuro* 2023 Sep PMID: 37775311
45. Wu SJ, Sevier E, Dwivedi D, et al. *Neuron* 2023 Sep PMID: 37390821
46. Chatzifrangkeskou M, Kouis P, Skourides PA. *J Cell Biol* 2023 Nov PMID: 37851005
47. Martins-Marques T, Witschas K, Ribeiro I, et al. *Open Biol* 2023 Nov PMID: 37907090
48. Subramanian M, Mills WT, Paranjpe MD, et al. *iScience* 2023 Dec PMID: 38235335
49. Christofidou, E;Tomazou, M;Voutouri, C, et al. *bioRxiv* 2023
50. Hayashi M, Kawasaki T, Nakatsukasa K. *FEMS Yeast Res* 2023 Dec PMID: 38130235
51. Do JL, Pedroarena-Leal N, Louie M, et al. *Invest Ophthalmol Vis Sci* 2023 Dec PMID: 38117244
52. Große-Segerath L, Follert P, Behnke K et al. *Nat Commun* 2024 Feb PMID: 38316785

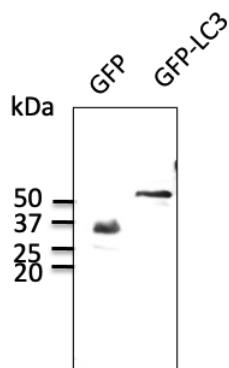


Anti-GFP Ab at 1/1,000 dilution; HEK293 transfected cell lysates at 50 µg per lane; rabbit polyclonal to goat IgG (HRP) at 1/10,000 dilution;

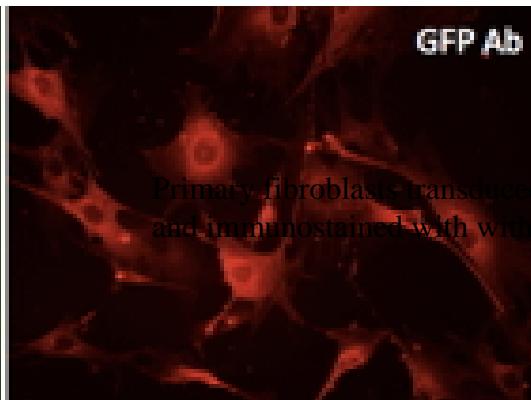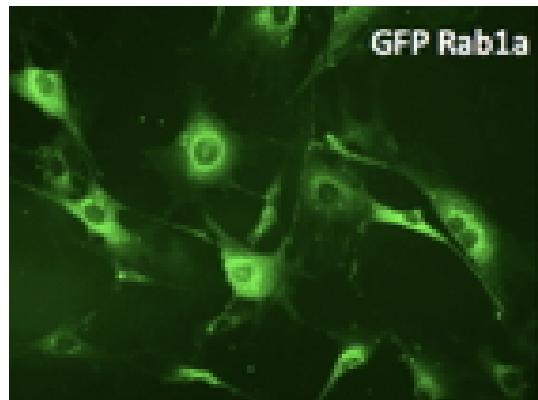

Primary fibroblasts transduced with Ad GFP-Rab1a and immunostained with GFP Ab;

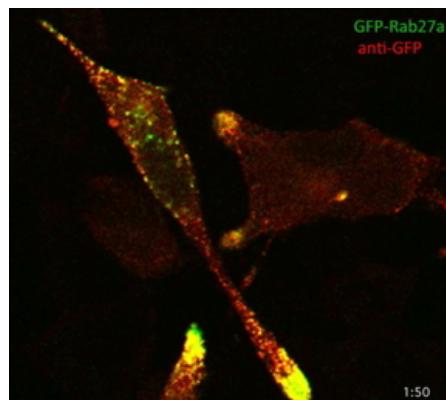

AtT-20 transfected with GFP-Rab27a and immunostained with GFP Ab;

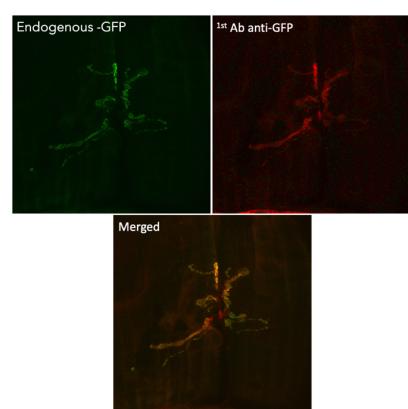

Immunofluorescence in Drosophila larvae NMJ muscle 6/7 expressing GluRIIA-GFP in neurons (GluRIIA is a post-synaptic protein) using 1st Ab anti-GFP at 1/1,000 and 2nd Ab anti-goat IgY conjugated to DyLight®550 at 1/500;

### **SICGEN's Proprietary Immunogen Policy**

In order to produce high specific antibodies SICGEN has invested a lot of time and effort into selecting immunogen sequences. SICGEN has decided to protect this information by not publishing it on the website. However, these sequences are available on request.